



When to propose the Stepper Line

The Stepper Orthoses have been designed to accommodate a **swollen and painful foot after surgery**. They provide support where needed while also offering **comfort and reassurance**. They protect injuries from accidental impacts, lighten the load, and instill a **pleasant sense of security** for walking stability.

9000P00202INT01.0424

stepper

M23001 • Stepper SURG post-op orthosis for the foot

| | | | | | | |
|-----------|---------|---------|-------|---------|---------|---------|
| Code | M23001 | | | | | |
| Size | XXS | XS | S | M | L | XL |
| Shoe size | 34/36,5 | 37/38,5 | 39/41 | 41,5/43 | 43,5/45 | 45,5/47 |

Supplied individually- ambidextrous

Indications

- Post-op for stabilization and reduce load from the foot (e.g. foot osteotomies, hallux valgus,...)

M23002 • Stepper DIGIT post-op orthosis to reduce load on the forefoot

| | | | | | | |
|-----------|---------|-----------|---------|-------|---------|--|
| Code | M23002 | | | | | |
| Size | XS | S | M | L | XL | |
| Shoe size | 33,5/37 | 37,5/39,5 | 40/41,5 | 42/44 | 44,5/47 | |

Supplied individually- ambidextrous

Indications

- Post-op orthosis to reduce load on the forefoot

M23003 • Stepper CAL Post-op orthosis to reduce load on the hindfoot

| | | | | | | |
|-----------|---------|-----------|---------|-------|---------|--|
| Code | M23003 | | | | | |
| Size | XS | S | M | L | XL | |
| Shoe size | 33,5/37 | 37,5/39,5 | 40/41,5 | 42/44 | 44,5/47 | |

Supplied individually- ambidextrous

Indications

- Post-operative recovery after both soft tissue and bony heel surgery
- Selective pressure relief after fractures
- Restoration of hindfoot ulcers



Stepper SURG



Stepper DIGIT



Stepper CAL

RADIO LUCENT


RIGID SOLES

SINGLE PATIENT MULTIPLE USE

UNISEX

stepper

Post-op orthoses for foot

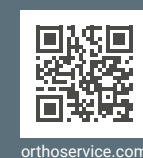

orthoservice.com



Headquarter: ORTHOSERVICE AG
 Via Milano 7 · CH-6830 Chiasso (TI) · Switzerland
 Tel. 0041 (0) 91 822 00 88 · Fax 0041 (0) 91 822 00 89
 info@orthoservice.com · www.orthoservice.com

ORTHOSERVICE DEUTSCHLAND GmbH
 Flugstraße 8 · D-76532 Baden-Baden · Deutschland
 Tel. 0049 (0) 7221 991 39 11 · Fax 0049 (0) 7221 991 39 13
 info@orthoservice.de · www.orthoservice.de

ORTHOSERVICE POLSKA Sp. z o.o.
 ul. Warszawska 416a · 42-209 Częstochowa · Polska
 Tel. +48 (0) 34 340 13 10 · (NIP): 9492246785
 info@orthoservice.pl · www.orthoservice.pl

RO+TEN S.r.l.
 Sede legale: Via Marco De Marchi 7 · I-20121 Milano (MI) · Italia
 Sede operativa e amministrativa:
 Via Comasina 111 · I-20843 Verano Brianza (MB) · Italia
 Tel. 0039 039 601 40 94 · Fax 0039 039 601 42 34
 info@roplusten.com · www.roplusten.com

ORTHOSERVICE RO+TEN
 Take care feel better

ORTHOSERVICE RO+TEN
 Take care feel better

Why choose these orthoses for post-op

- 1 Maximum functionality** thanks to the ultra-rigid sole
- 2 Maximum protection** thanks to the squared-off front shape to protect the toes
- 3 Maximum comfort** thanks to the lightweight and breathable upper and the wide opening for fit
- 4 Unisex**, which the warehouse thanks for



UNISEX



FULL EVA
INSOLE

TOTAL
OPENING



stepper SURG

- No height compensation necessary on the contralateral limb, low-profile
- Reduction of pressure on the heel and forefoot thanks to the rigid sole
- Boat-shaped profile that facilitates foot roll-off

UPPER

- Lightweight and breathable
- Rear counter for better hindfoot stabilization
- Fully openable for easy fitting
- Velcro® strap on the instep with reinforcement ring; with protective padding
- Symmetric closure on the forefoot, without buckles

RIGID SOLE

- Stiff non-slip sole with reinforcement bar
- Square-shaped front for toe protection
- Removable EVA insole to be optionally adapted or replaced

stepper DIGIT

- Wedge at 15° to reduce load on the forefoot



stepper CAL

- Selective unloading of the pressure on the heel



Exclusive
design by
Orthoservice RO+Ten